

















Index

# 01. Presentation 02. Fairs 04. Granulometry 05. Sectors and application examples 06. Sieves 08. Sieves and accessories 09. Ultrasonic bath: NEPTUNO 10. Mesh for laboratory sieves 11. Testing laboratory accredited by ENAC 12. Filters 14. Industrial mesh 15. Laboratory equipment 16. Digital electromagnetic sieve shaker: IRIS 20. Tapping sieve shaker: SALOMÓN 21. Digital air jet sieve shaker: EOLO 22. Laboratory grinding mill: ARES 24. Laboratory V mixer: VENUS 25. Industrial equipment 26. Circular sieve shaker: ZEUS 28. Circular side engines sieve shaker: HELIOS 29. Circular sieve shaker: VULCANO 30. Mobile sieve shaker: CENTAURO 31. Centrifugal sieve shaker: MINERVA 32. Industrial V-Mixer: ATENEA 33. Accessories and spare parts

FILTRA reserves the right to change the dimensions and technical characteristics of this catalog without notice.

# Presentation



FILTRA VIBRACION, S.L. has been proud to celebrate its **25 years** in the market developing and manufacturing laboratory and industrial process equipment for the sieving of solids and liquids and for the mixing and milling of solids. Another production line is the industrial filters made-to-measure according to customer specifications.

FILTRA offers high quality products at very competitive prices. All products strictly comply with national and international manufacturing standards, subjected to rigorous quality control and being delivered with an EC Declaration of compliance. The FILTRA quality management system for sieve manufacture is certified by TÜV NORD.

Following our policy of Quality and Innovation, FILTRA continues to investigate improvements and the development of new products.







Manufacture of stainless steel mesh or perforated plate test sieves of diameter between 60 and 600 mm.



Laboratory accredited by ENAC for the emission of reports of wire mesh and perforated plate test sieves according to the following standards: UNE 7050-3, ISO 3310-1, ASTM E11, ISO 9044, UNE 7050-4, ISO 3310-2, ASTM E323.



FILTRA is a member of the Association of Scientific Equipment Manufacturers LABMAS SPAIN, which collaborates with other international industry associations.



# Fairs



We distribute in more than 50 countries

We are leaders in the Spanish market, and currently sell in over 50 countries to multiple sectors -pharmacy, food, chemical industry, construction, recycling, agriculture, etc. - and applications.

Our hope is to continue growing and share our enthusiasm and satisfaction wiht our partners and customers.



Expoquimia Fair. Barcelona



Achema Fair. Munich









Powtech Fair. Nuremberg



Arablab Fair. Dubai



Exposólidos Fair. Barcelona



# Granulometry

Sectors and examples of application

Granulometric classification or granulometry is defined as the measurement and grading which is performed on granules or particles of a solid material for analytical purposes, both on its origin and its mechanical properties, and the calculation of the proportion of the corresponding elements for each size on a granulometric

### What is the significance of granulometry?

Knowledge of the granule dimensions and the distribution according to size of the particles dispersed in solids is extremely important for research and development, production and quality control applications.

Based on the granulometric distribution, some of the product properties will vary, such as:

- Flavour
- Absorption
- Conductivity
- Miscibility
- Agglomeration due to cohesion/adhesion forces

- Mechanical properties
- Fracture properties
- Resistance to wear
- Surface reaction
- Insulation properties
- Filtration properties



The FILTRA sieves shakers can be used to classify and analyse a wide variety of products, both liquid and solid.

This equipment can be used in areas of research and development, quality control, and production process supervision.

### The main sectors for application are:

**Construction** 

derivatives

**Cosmetics and** 

Plaster

<u>hygiene</u>

Soaps

**Mining** 

Detergents

**Pigments** 

Surfactants

Sand

Clay Cement and

Filling materials

Flours Milk products

Tiger nut milk Musts Cereals Spices

Pulses Chocolate Nuts Meat products Sugar

Coffee Mineral salts Fruit compots Ice-creams Juices

**Ceramics** Stains Enamels **Pigments** Clay

**Agriculture** Fertilizers Seed Grain Herbs Flowers

Chemical and **Petrochemical** Glues/adhesives **Paints** Toner Plastic and polyamides, etc Masterbach Pet

Recycling Paper pulp Waste water Glass Slurry Waste oils Plastics Copper

Metals Metal oxides Kaolin **Explosives** Minerals Diamonds

> **Pharmaceuticals** Medicines and drugs



# Sieves

FILTRA VIBRACION, S.L manufactures a wide range of laboratory sieves with diameters from 60 to 600 mm, with top quality AISI 304 or AISI 316 stainless steel for the frame, AISI 304 for the perforated plates, and AISI 316 for the metal meshes, providing excellent durability and high resistance to corrosion.



They are all manufactured in strict compliance with UNE, ISO, ASTM, AFNOR, BS, etc. national and international regulations. The FILTRA Quality Management System for the manufacture of these sieves is certified by TÜV NORD, in accordance with UNE-EN ISO 9001:2000 regulations.

They are manufactured in two pieces, allowing the mesh to be changed when it has deteriorated or does not meet the regulations. The interior sieve surface is totally flat to prevent the accumulation of residue and dirt. Cleaning is very fast and easy.



Manufacture of stainless steel mesh or perforated plate test sieves of diameter between 60 and 600 mm.

Our laboratory sieves have a rubber seal included in order to ensure they fit tightly. In addition, they can be fitted together and fitted to those from other manufacturers perfectly.

Marking is done by laser etching on the frame, giving each sieve an individual identification number, which allows for subsequent tracing.




### **QUALITY**

Before being supplied, each sieve is subjected to optical calibration and provided with a manufacturer's certificate and documentation regarding maintenance.

We also have lids and bottoms in stainless steel in all diameters for which we manufacture the sieves.



### Table of the various diameters and heights of FILTRA sieves

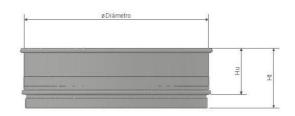

D	Hu	Ht	P	v
60	25	30	60	5
60	75	85	120	15
75	25	40	85	25
100	20	35	100	50
100	50	60	150	75
125	20	34	165	85
150	30	43	200	100
150	50	60	250	120
200	25	50	450	70
200	50	70	500	140
200	100	120	800	250
200	200	220	1000	500
203	25	50	450	70
203	50	70	500	140
250	70	90	900	280
300	80	95	1250	560
305	50	70	1250	560
315	80	95	1300	560
350	80	100	1500	600
400	65	85	1700	600
400	100	115	2000	650
450	100	115	2200	750
500	100	115	3000	850
600	110	125	3600	1000

D	DIAMETER (mm)
Hu	USEFUL HEIGHT (mm)
Ht	TOTAL HEIGHT (mm)
P THEORICAL WEIGHT (g)	
	(Depending on Mesh)
٧	RECOMMENDED VOLUME
	(with mesh of 1 mm cm <sup>3</sup> )

### **ADVANTAGES:**

- Manufacture according to national and international regulations.
- High quality. High resistance to corrosion. Excellent durability.
- Manufactured in two pieces: allows the mesh to be changed.
- Totally flat surface for easy cleaning. Product residues do not accumulate.
- Individual identification by laser, allows the sieve to be traced.
- The sieves can be fitted together and fitted to those from other manufacturers perfectly.
- Maximum stability and sealed fit for sieves stacks.



 $\epsilon$ 

# **Ultrasonic Bath**

FILTRA offers a wide range of sieves in a variety of sizes, shapes and materials, together with complementary equipment.



Sieve with stainless steel frame and mesh.



Stainless steel sieve with square-holed perforated plate.



Stainless steel framed sieve with nylon mesh.



Stainless steel sieve with oval-holed perforated plate (for cereals, ISO 5223).



Stainless steel sieve with round-holed perforated plate.



Special stainless steel sieve for EOLO-FTLBA (Air Jet) sieve shaker, with stainless steel or nylon mesh.



Special test sieves: conical, double height, half height, etc.



Grid Sieves made of paint carbon steel.



Stainless steel lids and receivers in all sizes.



Sieve shelves.



Cleaning brushes for sieves.



Food balls used to open the mesh of sieve.

A NEPTUNO ultrasonic bath is the perfect accessory to maintain and clean laboratory sieves and other pieces of equipment.

Ultrasonic cleaning is currently the most modern and effective cleaning system: the alternate high and low pressure waves produce a frequency which, through the cavitation phenomenon, leads to the production of millions of microscopic bubbles which impact on the object surface, loosening dirt in minimum time.


Ultrasonic baths have low electrical energy consumption and save time when compared to brushing out or other traditional methods. The areas which are difficult to access become clean and without traces of dirt, while the object surface is not damaged, as can happen with other methods.

The NEPTUNO ultrasonic baths are made from stainless steel (AISI 304). We have various models according to capacity: 0.7/1/1.5/3/6/10/14/28/35/50 L.







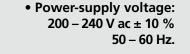

**GENERAL EQUIPMENT** 

**CHARACTERISTICS:** 

- Current and power surge protection.
- Digital control panel. It allows monitoring of the receptacle temperature, modification of the operating temperature, ultrasonic action time and the switching on and activation of the cleaning equipment.
  - High performance steelaluminium piezoelectric transducers (the number of transducers will depend on the power of the equipment).
    - Working temperature range: from 0 to 80 °C.



	Frequency Ultrasound (KHz)	Power Ultrasound (W)	Power Heating (W)	Dimensions Internal (mm) W x L x H	Weight (Kg)	Drainage
COMBAÑO 0,7L	40	50	100	85 x 150 x 60	1,8	NO
COMBAÑO 1L	40	50	100	120 x 130 x 65	2,1	1/2".
COMBAÑO 1,5L	40	50	100	130 x 140 x 100	2,3	1/2".
COMBAÑO 3L	40	100	100	130 x 230 x 100	3,7	1/2".
COMBAÑO 6L	40	200	150	150 x 300 x 150	5,7	1/2".
COMBAÑO 10L	40	300	200	235 x 295 x 150	9	1/2".
COMBAÑO 14L	40	500	300	295 x 325 x150	14	1/2".
COMBAÑO 28L	40	600	750	295 x 495 x 200	18	1/2".
COMBAÑO 35L	40	700	1000	295 x 395 x 300	25	1/2".
COMBAÑO 50L	40	800	1000	300 x 600 x 300	40	1/2".

OB CENTRAL CE

# Mesh / Perforated Plate for Laboratory Sieves

Table of normalized openings in metal mesh and plate for the manufacture of laboratory sieves, according to various international Regulations

### **METALLIC MESH**

UNE 7050-3 ISO 3310-1	AST	TM E11
W*	W*	N°
125	125	5 in.
112		1
106	106	4,24 in.
100	100	4 in.
90	90	31/2 in.
80		İ
75	75	3 in.
71		İ
63	63	21/2 in.
56		İ
53	53	2,12 in.
50	50	2 in.
45	45	13/4 in.
40		İ
37,5	37,5	11/2 in.
35,5		
31,5	31,5	11/4 in.
28		İ
26,5	26,5	1,06 in
25	25	1 in.
22,4	22,4	7/8 in.
20		
19	19	3/4 in.
18		1
16	16	5/8 in.
14		
13,2	13,2	0,530 in.
12,5	12,5	½ in.
11,2	11,2	7/16 in.
10		
9,5	9,5	3/8 in.
9		
8	8	5/16 in.
7,1		
6,7	6,7	0.265 in.
6,3	6,3	1/4 in.
5,6	5,6	31/2
5		
4,75	4,75	4
4,5		
4	4	5
3,55		
3,35	3,35	6
3,15		1
2,8	2,8	7
2,5		
2,36	2,36	8
2,24		
2	2	10
1,8		-
1,7	1,7	12
1,6		1
1,4	1,4	14
1,25		1
1,18	1,18	16
1,12		1
1	1	18

UNE 7050-3 ISO 3310-1	ASTM E11	
W*	W*	N°
900		
850	850	20
800		
710	710	25
630		
600	600	30
560		
500	500	35
450		
425	425	40
400		
355	355	45
315		
300	300	50
280		
250	250	60
224		
212	212	70
200		
180	180	80
160		
150	150	100
140		
125	125	120
112		
106	106	140
100		
90	90	170
80		
75	75	200
71		
63	63	230
56		
53	53	270
50		
45	45	325
40		
38	38	400
36		
32	32	450
25	25	500
20	20	635

### PERFORATED PLATE

mm (milimeters)

UNE 7050-4 ISO 3310-2	ASTM E323		
W*	W*	N°	
● ■ 125	● ■ 125	5 in.	
● <b>■</b> 112		† 1	
● ■ 106	● ■ 106	41/4 in.	
● ■ 100	● ■ 100	4 in.	
● ■ 90	● ■ 90	31/2 in.	
● ■ 80	1 -		
● ■ 75	● ■ 75	3 in.	
<b>● ■</b> 71	1		
● ■ 63	● ■ 63	21/2 in.	
● ■ 56			
● ■ 53	● ■ 53	21/8 in	
● ■ 50	● ■ 50	2 in.	
● ■ 45	● ■ 45	13/4 in.	
● ■ 40			
● ■ 37,5	● ■ 37,5	11/2 in.	
● ■ 35,5	İ		
● ■ 31,5	● ■ 31,5	11/4 in.	
● ■ 28	1	į i	
● ■ 26,5	● ■ 26,5	11/16 in.	
<b>●</b> ■ 25	● ■ 25	1 in.	
● ■ 22,4	● ■ 22,4	7/8 in.	
● ■ 20			
● ■ 19	● ■ 19	3/4 in.	
● ■ 18			
● ■ 16	● ■ 16	5/8 in.	
● ■ 14			
● ■ 13,2	● ■ 13,2	17/32 in.	
● ■ 12,5	● ■ 12,5	1/2 in.	
● ■ 11,2	● ■ 11,2	7/16 in.	
● ■ 10			
● ■ 9,5	● ■ 9,5	3/8 in.	
● ■ 9			
● ■ 8	● ■ 8	5/16 in.	
● ■ 7,1			
● ■ 6,7	● ■ 6,7	17/64 in.	
● ■ 6,3	● ■ 6,3	1/4 in.	
● ■ 5,6	● ■ 5,6	7/32 in.	
● ■ 5	ļ		
● ■ 4,75	● ■ 4,75	3/16 in.	
● ■ 4,5			
<b>●</b> ■ 4	● ■ 4	5/32 in.	
● 3,55	1 2 2 2 5	4/0:	
● 3,35	● 3,35	1/8 in.	
● 3,15	100	7/64 :	
● 2,8	● 2,8	7/64 in.	
● 2,5 ● 2,36	<b>■</b> 2.36	3/62 in.	
● 2,36 ● 2,24	● 2,36	3/02 III.	
• 2	•2	0,078 in.	
● 1,8	•	0,070 III.	
● 1,7	● 1,7	0.066 in.	
● 1,6	1		
● 1,4	● 1,4	0,055 in.	
● 1,25	1		
● 1,18	● 1,18	0,045 in.	
● 1,12	†		
<b>●</b> 1	<b>●</b> 1	0,039 in.	
	1 -		

# ENAC (Spanish National Accreditation Agency) Test Laboratory

In 2002 FILTRA VIBRACION, S.L set up the first ENAC accredited laboratory in Spain (Accreditation no. 310/LE683) for the carrying out of Trials for metal mesh and plates in accordance with the following regulations:

Mesh	Plates
UNE 7050-3	UNE 7050-4
ISO 3310-1	ISO 3310-2
ASTM E11	ASTM E323
ISO 9044	

The laboratory technical personnel, using artificial vision equipment and following methods established by the UNE-EN ISO/IEC 17025:2005 regulation governing our quality systems, carry out necessary measurements and calculations in order to issue a trial report which, together with other results, determines whether the mesh or perforated plate of the examined sieve conforms to the reference regulation.

ENAC accreditation is recognized internationally through MLA (Multilateral Agreement) or MRA (Mutual Recognition Agreement).







\*W = Mesh light (opening size)

# **Filters**

### **METALLIC MESH FILTERS**

At FILTRA VIBRACIÓN, S.L. we specialize in the manufacture of all types of metal mesh industrial filters. We have a team with over 20 years' experience in this field, developing, together with our customers, solutions in filter consumables for multiple industrial sectors such as chemicals, petrochemicals, pharmaceutical, plastics, resins, thermal oils, etc.

We manufacture a multitude of industrial filters in AISI 304 or AISI 306 stainless steel: cylindrical filters, DN22 to DN 600Y-type; DN25 to DN600 filter baskets, pleated filters with one or various mesh layers; conical and truncated conical filters; interchangeable plugs; mesh disks; special elements.



We develop and manufacture made-to-measure products according to the client's plans or specifications, from a single unit to whole ranges.














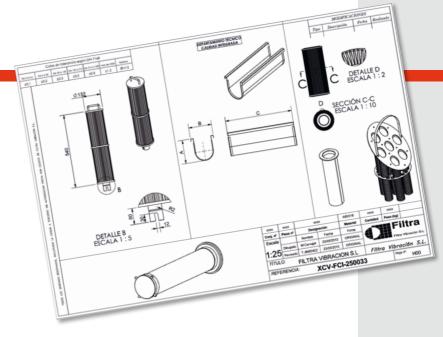

### **METALLIC MESH SieveS**

We manufacture all types of custom metal sieves with a single mesh or multiple or with reinforcement structure where necessary.





### **NON-METALLIC FILTERS**

We supply synthetic filters made of polypropylene, polyester or polyamide: sleeves and bag filters, mini weaves, filter discs, cartridges, mesh fabrics for fluidization and nylon mesh bulk meters.









# **Industrial Mesh**

### **METALLIC MESH**

We have a large stock of stainless steel mesh AISI 304 and 316 acc/ISO 9044 to supply in rolls, panels or bulk meters: plain fabric, twill, wave, electro, sintered or Reps.




### Table of standardized opening mesh sizes for production of industry sieves.

Strong Group


Lightweight Group			
Opening Ø Wire % size (mm) (mm) Surface		Useful	
1,440	0,300	68,80	
1,126	0,280	66,50	
1,150	0,240	68,50	
1,030	0,230	66,50	
0,880	0,230	62,40	
0,726	0,200	62,00	
0,594	0,200	55,20	
0,534	0,160	58,20	
0,457	0,160	55,50	
0,416	0,140	57,10	
0,385	0,120	56,50	
0,343	0,120	53,80	
0,307	0,120	52,50	
0,277	0,120	49,70	
0,270	0,100	53,10	
0,247	0,100	51,80	
0,219	0,090	50,60	
0,198	0,080	51,80	
0,188	0,090	45,60	
0,173	0,080	46,30	
0,161	0,070	48,30	
0,150	0,063	49,30	
0,128	0,070	49,30	
0,125	0,060	45,50	
0,124	0,050	50,10	
0,104	0,050	45,20	
0,089	0,050	40,90	
0,080	0,050	46,30	
0,075	0,036	36,20	
0,063	0,063	34,80	
0,056	0,036	35,30	
0,049	0,030	34,80	
0,042	0,036	28,00	
0,040	0,030	33,20	
0,025	0,025	24,90	

Opening size (mm)	Ø Wire (mm)	% Useful surface
19,716	2,5	78,80
11,880	2	73,40
9,5	1,6	73,40
7,96	1,3	73,70
5,84	1,1	70,30
4,56	1	67,00
4	1	64,30
3,73	0,9	64,80
3,17	0,8	63,60
2,77	0,7	63,50
2,39	0,7	59,40
2,18	0,6	61,40
1,81	0,5	60,90
1,48	0,5	55,50
1,3	0,44	54,30
1,12	0,42	52,10
0,99	0,40	50,60
0,75	0,36	45,50
0,606	0,32	43,20
0,514	0,28	41,20
0,454	0,24	41,90
0,36	0,20	41,30
0,283	0,18	36,40
0,237	0,16	36,40
0,207	0,14	36,40
0,189	0,12	37,00
0,168	0,11	36,40
0,153	0,10	36,10
0,141	0,09	37,00
0,125	0,08	35,20
0,111	0,063	40,00
0,104	0,050	45,63
0,088	0,050	40,75
0,076	0,045	36,55
0,061	0,045	30,25
0,057	0,035	38,30
0,049	0,030	31,00
0,039	0,030	33,70
0,032	0,025	31,00
0,025	0,025	25,00
0,020	0,020	25,00


Superstrong Group				
Opening size (mm)	Ø Wire (mm)	% Useful surface		
25	3	79,50		
19,72	2,5	78,80		
16,51	2	79,40		
10	2,2	67,00		
8	2	73,00		
6,3	1,6	64,00		
5	2	51,00		
5	1,2	65,00		
4	1	64,00		
3	1	56,00		
2,5	1	51,00		
2	1	44,00		
1,6	0,70	48,00		
1,5	1	36,00		
1,23	0,5	51,00		
1,18	0,55	47,00		
1,08	0,65	39,00		
0,9	0,36	51,00		
0,82	0,44	42,00		
0,71	0,40	41,00		
0,52	0,40	32,00		
0,40	0,22	42,00		
0,30	0,20	36,00		

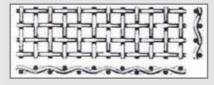


Plain fabric square

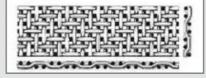

Plain fabric twill

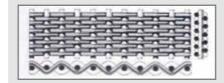

Fabric Reps together

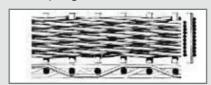

Fabric Reps twill together



# Electromagnetic Sieve Shaker

# IRIS

# IRIS

The FILTRA electromagnetic sieve shakers are used in research and development, and in quality control for production processes, both for raw materials and intermediate and finished products. With them it is possible to carry out granulometric testing for particle separation, both for dry and wet conditions.

The IRIS model sieve shakers are highly effective, thanks to the three-dimensional particle motion, created by an electromagnet. With this 3D effect the particles bounce over the whole sieve area, forcing them to pass through to the lower sieve more quickly. Thus a faster, more effective sieving is achieved through the whole sieve column.

There are four models: FTS-200, FTL-200, FTL-300 and FTL-400, which cover sieves from 60 mm Ø to 400 mm Ø for samples from 25 gr to 25 kg of product. The particle size to be analysed can range from 25 $\mu$  to125 mm. Each user can choose the equipment which is best suited to their needs.

### **Examples of application:**

Foodstuffs, earth, coal, construction materials, chemical and pharmaceutical products.



# Filtre

### IRIS FTL-0200

Digital electromatic sieve shaker for sieves from Ø 60 a 203 mm. Its three-dimensional movement throughout the sieve column makes it very effective when carrying out granulometric test. It includes digital power, timing and intermittence controls for sieving. It can be supplied with the standard anchoring system or with the "Easy-Press" rapid system. It permits 8 sieves of 50 mm in height, or 16 sieves of 25 mm in height.

### **ANCHORING SYSTEMS:**

 Standard: Threaded rods, bakelite knobs and methacrylate lid.



 "Easy-Press": smooth rods, pressurized tightening knobs which easily slide along the rods, and a methacrylate lid.








All models include digital power control (from 1 to 9), sieving time (from 0 to 99 min.) and the sieving mode; continuous or intermittent.

All the sieve shakers are equipped with an anti-vibration system which, together with the anchoring system, ensures complete stability during the sieving.

Depending on the model, it is possible to choose between the standard fixation system, and the "Easy-Press" rapid fixation system. The methacrylate lid makes observing the product during the sieving process easy.

All equipment is subjected to strict quality control, which guarantees years and years of problem-free use.

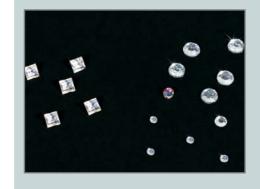

### **IRIS FTS-0200**

Small, compact and light, it is the economic option from the IRIS range of electromagnetic sieves. Allows sieves from Ø 60 to 203 mm Up to 7 sieves of 50 mm in height, or 14 sieves of 25 mm in height. It includes digital power, timing and intermittence controls for sieving. It is supplied with the standard fixation system.



### **ADVANTAGES:**

- Three-dimensional motion.
   Highly effective: short sieving times
- Sieving for dry and wet test
- Ergonomic design
- Digital adjustment of sieving parameters
- Easy to operate
- Integrated interface
- Anti-vibration system
- Silent
- Maintenance free

16 CE

## **Examples of application:**

Minerals, fertilizers, plastics, metals, seeds, flowers, herbs, grain, woodchip, fibres, etc.





4		13	200	(3)
10	7	<b>M</b>		120
10		Ŧ.	1	John
2	N.	7		
50	73			
			No.	M

# FTL-0400



Designed for high sample volumes, it can admit 5 sieves of 80 mm in height, or 8 sieves of 50 mm in height. The diameters range from 250 to 305 mm. The power, time of sieving and intermittence are all digitally adjusted.

It can be supplied with the standard fixation system or with the "Easy-Press" rapid system, according to your choice.

### **ACCESSORIES:**

- Methacrylate lid for dry test.
- Methacrylate lid with stainless steel diffusing ball for wet test.



- Stainless steel receiver for dry test.
- Stainless steel receiver with outlet for wet test (including hose and flange).
- Stainless steel lid.
- "Easy-Press"rapid anchoring system (smooth rods and pressurized tightening knobs).
- Feed balls to ease the passage of the product in the sieve.
- Cleaning brushes.
- Freeing tool.



IRIS FTS-0200 FTL-0200 FTL-0300 Sieve Shaker Type Electromagnetic Electromagnetic Electromagnetic Electromagnetic Dry test Yes Yes Yes Yes Wet test Yes Yes Yes **Levels of Power** 9 Size of Particle 20 μ a 125 mm 20 μ a 125 mm 20 μ a 125 mm 20 μ a 125 mm 5 x Ø300x80 7 x Ø200x50 8 x Ø200x50 4 x Ø400x100 Number of Sieves (mm) 14 x Ø200x25 16 x Ø200x25 8 x Ø300x50 315 Max. Diameter of Sieve (mm) 203 203 400 Timer (min) 1 a 99 1 a 99 1 a 99 1 a 99 220/240 V - 50/60 Hz. Tension 220/240 V - 50/60 Hz. 220/240 V - 50/60 Hz. 220/240 V - 50/60 Hz. Intensity (A) 0,8 0,8 0,8 0,8 Power (kW) 0,4 Weight (Kg) 27 40 45 48 305 x 630 254 x 800 254 x 800 555 x 630 x 750 Dimensions L x W x H (mm)

IRIS FTL-0400

controls for sieving.

system.

the standard anchoring system

# Tapping Sieve Shaker

# SALOMON

# Digital Air Jet Sieve Shaker

## EOLO

### **FT-VG-200**

The Sieve Shaker Salomon performs the sieving process reproducing the movements of the hand sieving, the combination of oscillating movements and tapping on the top part facilitates the passage of the particles of the sample, avoiding the clogging of sieves.

The equipment is fully provided with acoustic insulation to minimize the degree of noise during sieving.



### Inclues:

- Sound Enclosure Cabinet provided with acoustic insulation and protection sensor to stop the machine when the door is open for Safety.
- Metal support for sieves with adjustable height.
- Top metal lid with rubber stopper.
- Tapping arm situated in the upper part of equipment, driven by an axis which acts jointly with the oscillating movement of sieving and allows the beating at the top of the sieves.
- Driving set composed of a motor of 1.5 kW at 220 V monophase and the corresponding pulleys toothed belt transmission.

TECHNICAL CHARACTERISTICS		
Number of Sieves (mm) 6 x Ø 200 x 50		
Rotation speed (r.p.m.)	275	
Tapping speed (taps/min)	144	
Dimensions L x W x H (mm)	810 x 540 x 750	
Electrical connection (V)	220 Monophase	
Power (W)	750	
Weight (Kg)	178	

### **Examples of application:**

Sands, chemicals and pharmaceuticals products, seeds, shavings, plastics, wood, etc.







# DRIVING DIGITAL MODULE INCLUDING:

- Switch on / off process.
- Switch "mushroom" emergency stop.
- Digital programmer sieving time.

### **FTLBA**

The EOLO digital air jet sieve shaker has been specially designed for dry granulometric testing, for fine particles from 5 microns in size (from 0.005 to 4 mm).

It offers fast, effective sieving thanks to the flow of the product obtained forced by a current of air which forces the particles through the sieve. This effect is obtained using an industrial vacuum which maintains constant depression.


It features a digital panel with sieving time regulation (1 to 99 minutes) and Vacuum Meter, which incorporates a depression regulator valve in order to select the necessary depression for each test.

The sieve shaker includes a methacrylate lid which allows the product behaviour to be observed, and a nylon mallet to loosen the residual product which is stuck to the lid due to the effects of static electricity.

The use of FILTRA quality Air Jet sieves of Ø200 mm (not included) is recommended.



TECHNICAL CH	ARACTERISTICS
Sieve Shaker Type	Air jet
Dry test	YES
Wet test	NO
Size of Particle	5 μm a 4 mm
Number of Sieves	1
Diameter of Sieve (mm)	200
Timer (min)	1 a 99
Tension	220/240 V - 50/60 Hz.
Power (kW)	0,19
Máximum Depression (Kpa)	20
Weight (Kg)	20
Dimensions ØxH (mm)	360 x 275

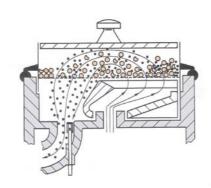

Current of air which forces the particles through the sieve.

### **Examples of application:**

Food powders, fine sand, construction materials, chemical and pharmaceutical products, minerals, fertilizers, powdered plastics, powdered metals, fibres, etc.



### **ACCESSORIES:**

- Methacrylate lid.
- Nylon mallet.
- NILFISK GM-80 vacuum.
- FILTRA Ø 200 mm Air Jet sieves.
- Cleaning brushes.
- Freeing tool.
- Ultrasound bath for sieves cleaning.



### **ADVANTAGES:**

- Highly effective in the sieving of fine particles: very short sieving times.
- Dry sieving.
- Ergonomic design.
- Digital adjustment of sieving parameters.
- Easy to operate.
- Integrated interface.
- Silent
- Maintenance free

20 Filtra

### FML-0100 / FML-2000

FILTRA manufacture two different models of laboratory Grinding Hammer Mills, in order to cover the different customer needs: ARES FML-0100 is a small and economical one for milling small quantities of product and ARES FML-2000 with the double of capacity for biggest quantities of product or biggest articles.

They are used to grind a wide variety of soft, semi-hard and hard materials, such us, cereals, spices, vegetables, plastics, resins, minerals, earth, stone, coal, etc. with a particle feed size of +/- 15 mm in FML-0100 and of +/- 20mm in FML-2000.

The feeding is supplied by a hopper located at the top and can be dosed manually using a chopper. The door has a quick locking with a safety micro-sensor which stops the rotation of the hammers when it's opened.

Both the front and loading hopper are manufactured from mirror-polished AISI 304 stainless steel. In the interior the 3 hammers are made from anti-wear AISI 410 stainless steel. The equipment includes an easily interchangeable sieve with round perforations.

In ARES FML-2000 model, plus the sieve, the blades and the toothed crown are also interchangeable, and can be easily removed for cleaning without using any key. This model incorporates a potentiometer to regulate the engine speed.

### **Examples of application:**

Grain food, cereals, seeds, spices, vegetables, stones, earth, minerals, resins, plastics, chemical products, construction materials, ceramics, pharmaceutical products, etc.







MODEL	Timer (min)	Tension	Power (Kw)	Speed (rpm)	Grinding chamber Dimensions (mm)	Loading hopper Volume (I)	Grinding chamber Volume (I)	Drawer Volume (I)	Weight (kg)	Number of hammers	Size of sieve holes (Ø mm)	Máximum Hardness (Mohs)	Maximum size of particle at the entrance (mm)	General dimensions (mm)
FML- 0100	0 - 99	220-240 V (50-60 Hz)	1,9	3000	Ø110x66	1,2	0,4	1	38	3	1,2,3,4,5	6-7	15	225 x600 x420
FML- 2000	0 - 99	220-240 V (50-60 Hz)	2,6	3000	Ø150x88	5,5	1,2	2	72	3	0.5, 0.8, 1, 1.5, 2, 2.5, 3,4,5,6	6-7	20	520 x700 x700

### **ARES FML-0100**













### **ADVANTAGES:**

- Easy cleaning thanks to the interchangeable toothed crown and sieves without the need of tools.
- Potentiometer to regulate motor velocity.
- Removable polished stainless steel loading hopper and collecting drawer, for ease of cleaning.
- Very fine end granulometry.

<sup>22</sup> Filtra

# Laboratory V-Mixer

# **VENUS**

### FTLMV-0,5 / 01 / 02 / 04 / 08

The VENUS V-shaped laboratory mixers are used in the precision homogenization of powdered or granulated products in the pharmaceutical, chemical and food industries. It can gently and evenly mix solids with solids, in any percentage, or solids with liquids (maximum 10%).

It has two loading inlets with silicon lids and one outlet with a dosing valve. The product introduced revolves through 360° via mechanical transmission. Its characteristic V-shape means axial currents are produced when the body is turned, and these separate and unite the material. This, along with the radial mixing action, produces an extremely homogeneous result (in the proportion of 1 part in 10.000) without the need for paddles or other devices.

The mix containers are manufactured from AISI-316 stainless steel with a polished mirror finish inside and outside. Its angle and corner-free design prevents the accumulation of product residue in the interior.

The mixer works with a single-phase motor and has a timer, stop/start switch and fuse box.

The useful capacity, in order to achieve the correct homogenization, is 50% of the total mixer volume. The FILTRA V mixers are available in 0'5, 1, 2, 4 and 8 litres (total capacity).



### **Examples of application:**

Food powder, chemicals and pharmaceuticals, metals, plastics, additives, etc.





### **ADVANTAGES:**

- Precision mixing. Maximum product homogenization.
- High quality materials and finishes.
- Free of interior angles or corners; there is no build up of product residue.
- Silent.
- Maintenance free

MODEL	Tin	ner	Tension	Mixer Power (Kw)	Weight (kg )	Useful Volume (I)	N° de vueltas deposito x min. (rpm)	Dimensions L x W x H (mm)		
FTLMV-0,5	Mechanical 60 min.		220-240 V 50-60 Hz	0,1	15	0,25	38	415 x 400 x 415		
FTLMV-01	Mechanical	60 min.	220-240 V 50-60 Hz	0,1	17	0,5	38	415 x 400 x 415		
FTLMV-02	Mechanical	60 min.	220-240 V 50-60 Hz	0,1	17	1	38	415 x 400 x 430		
FTLMV-04	Mechanical	60 min.	220-240 V 50-60 Hz	0,1	20	2	25	550 x 400 x 550		
FTLMV-08	Mechanical	60 min.	220-240 V 50-60 Hz	0,25	25	4	18	700 x 450 x 650		



Circular Sieve Shaker

**ZEUS** 

# ZEUS

### FTI-0550 / 0800 / 1200 / 1500

The ZEUS sieve shaker is suitable for large production volumes in industrial sieving. It can be used for separation processes for solids in a liquid, for the separation of different granulometries within the same product, in safety sieving for the elimination of foreign bodies, for the homogenization of particles, etc.

It permits a production volume of 100 to 25,000 kg/h, depending on the product, the mesh opening size and the feed flow to the machine. It is suitable for sieving from 36 microns to 25 mm.

Manufactured from AISI 304 stainless steel (or AISI 316 on request), it features an imbalanced-weight motor which produces a three-dimensional motion which makes the sieving of the product very effective. Adjustment of the imbalanced-weights allows the regulation of product behaviour on the mesh and the sieving speed.

The bottom has a side outlet through which the obtained product can be unloaded. The sieves can be ordered with or without product outlet. In the case of choosing a sieve with outlet, the sieving process can be automated: as the product enters the sieve shaker, it also leaves via the sieve opening.

The ZEUS model allows up to 4 sieves, with which up to 5 product separations are obtained. The possible diameters are: 550, 800, 1200 or 1500 mm. As an option it can be ordered with fast-closing clamps, lid, viewing windows in the lid and body, and wheels. If the sieve shaker is going to be used in areas classified as being potentially explosive environments then it is manufactured with ATEX materials, and FILTRA will supply the corresponding ATEX certificate. No special floor or foundation is necessary for its installation.





### **Examples of application:**

Chemical, pharmaceutical, agricultural, fishery, forestry, biological and food products, paints, sand, metals, wood, plastic, glass, recycling, waste water treatment, etc.






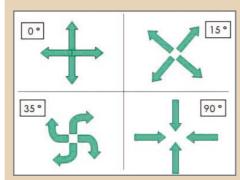

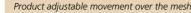



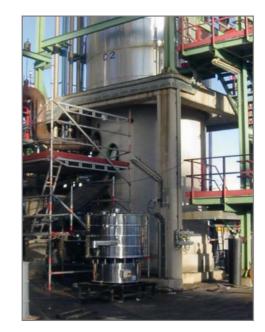









### **ADVANTAGES:**

- Sieve for large production volume.
- Highly robust. It can work for 24 hours.
- Highly versatile. Easy to operate.
- Sieves with or without product outlet (to choose).
- Can have wheels fitted
- Low maintenance.

MODEL	Sieve Shaker	Max. Number	Sieve	Tensi	on (V)	Intens	ity (A)	Powe	r (Kw)	Speed	(rpm)	Weight *	Dimensions
MIODEL	Туре	of sieves	Diameter	50 Hz	60 Hz	50 Hz	60 Hz	50 Hz	60 Hz	50 Hz	60 Hz	(Kg)	*ØxH mm
FTI-0550	Circular	4	550	400	460	1	1,02	0,42	0,5	1500	1800	95	710 * 755
FTI-0800	Circular	4	800	400	460	2,2	2	1,20	1,25	1500	1800	150	710 x 755
FTI-1200	Circular	4	1200	400	460	3	3,1	1,80	2	1500	1800	320	710 x 755
FTI-1500	Circular	4	1500	400	460	3	3,1	1,80	2	1500	1800	370	710 x 755

<sup>\*</sup> Size and weight for standard configuration: 1 level + lid

<sup>6</sup> Filtra (€

# HELIOS

# Circular Sieve Shaker

# **VULCANO**

### FTI-2M-0550 / 0800 / 1200

The HELIOS sieve shaker is especially recommended for control and safety sieving for both solids and liquids. It is very useful in confined spaces where the reception of material and the unloading of the sieved product are in close proximity. It permits a large product volume and a single sieving level.

They are often installed at the outlet of discharge chutes, grinding mills, mixers, etc. The inlet feed and product outlets are on the centre of the sieve shaker itself, it acting practically as a strainer. It is easy to install. It can be supplied on a wheeled support or alone, to be inserted within an installation.


It is made from stainless steel (AISI 304 or AISI 316) with two lateral motors, fitted on opposite sides, through which a three-dimensional motion in the horizontal and vertical planes is achieved, allowing excellent material flow.



If the sieve shaker is going to be used in areas classified as being potentially explosive environments then it is manufactured with ATEX materials, and FILTRA will supply the corresponding ATEX certificate.

It is available in three diameters: 550, 800 and 1200 mm.



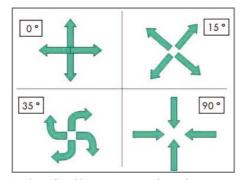

Product adjustable movement over the mesh.

### **Examples of application:**

Chemical, agricultural, food, pharmaceutical products, paints, resins, waste water, plastic, glass, etc.





### **ADVANTAGES:**

- Extremely robust and reliable. It can work 24 hours a day.
- It occupies little space.
- It allows for a large product volume.
- Central loading and unloading of material.
- Low maintenance.

### FTI-0300 / FTI-0400

The VULCANO Industrial Sieve shaker is ideal for small production runs of between 5 and 75 kg, depending on the type of product.

Manufactured from AISI 304 stainless steel (or AISI 316 on request), it features an imbalanced-weight motor which produces a three-dimensional motion which makes the sieving of the product very effective. Adjustment of the imbalanced-weights allows the regulation of product behaviour on the mesh and the sieving speed.

The bottom has a side outlet through which the obtained product can be unloaded. The sieves can be ordered with or without product outlet. In the case of choosing a sieve with outlet, the sieving process can be automated: as the product enters the sieve shaker it also leaves via the sieve opening.

This sieve shaker allows up to 5 sieves of 300 or 400 mm Ø. The standard fixation system is made up of threaded steel bars, bakelite knobs and an open (for continuous feed) or closed methacrylate lid. Optionally it can be ordered with a rapid closing system, metal lid and wheels for easy moving.

If the sieve shaker is going to be used in areas classified as being potentially explosive environments then it is manufactured with ATEX materials, and FILTRA will supply the corresponding ATEX certificate.







MODEL



60 Hz

50 Hz

60 Hz

50 Hz

60 Hz

### **Examples of application:**

Chemical, pharmaceutical and agricultural products, food stuffs, paint, sand, seeds, metals, wood, plastic, glass, etc.





### **ADVANTAGES:**

- Sieve shaker for small production volumes (pilot plants).
- Highly versatile.
   Easy to operate.
- Sieves with or without product outlet (to choose).

Dimensions
\*ØxH mm

- Can have wheels fitted
- Low maintenance.

Speed (rpm)

60 Hz

50 Hz

MODEL	Sieve Shaker	Max. Number	Sieve	Tensio	on (V)	Intens	ity (A)	Powe	r (Kw)	Speed	(rpm)	Weight *	Dimensions
	Туре	of sieves	Diameter	50 Hz	60 Hz	50 Hz	60 Hz	50 Hz	60 Hz	50 Hz	60 Hz	(Kg)	*ØxH mm
FTI- 2M-0550	2 motors	1	550	400	460	0,42	0,4	0,2	0,22	1500	1800	98	600 x 1000 x 615
FTI- 2M-0800	2 motors	1	800	400	460	0,42	0,4	0,2	0,22	1500	1800	150	800 x 1240 x 675
FTI- 2M-1200	2 motors	1	1200	400	460	0,63	0,66	0,36	0,4	1500	1800	210	1200 x 1705 x 815
* Size and w	eight for sta	ndard config	uration: 1 lev	el + cover									

-0800	ZIIIOLOIS	'	800	400	460	0,42	0,4	0,2	0,22	1500	1600	150	x 675		F11-0300	Circular	,	300	400	400	0,54	0,34	0,155	0,175	1500	1600	40	400 x 009	
TI- ·1200	2 motors	1	1200	400	460	0,63	0,66	0,36	0,4	1500	1800	210	1200 x 1705 x 815		FTI-0400	Circular	5	400	400	460	0,60	0,60	0,24	0,27	1500	1800	50	580 x 900	
e and v	veight for star	dard configu	ıration: 1 lev	rel + cover											* Size and we	eight for stan	dard configu	ıration: 1 leve	l + cover										

Filtra CE

### FTIS-0450 / FTIS-0600

The CENTAURO is a compact, light, mobile sieve shaker. It is frequently used to carry out on-site checking and safety sieving (mountains, fields, factory, etc.) as it can easily be taken from place to place. It can be used both for liquids and solids.

It is used for products such as: enamels, paints, ceramic powders and sands, in addition to applications in the chemicals, food and agricultural industries, etc.

It is made from AISI 304 stainless steel (it can also be made to order from AISI 316, and with materials suitable for ATEX explosive atmospheres, bearing the corresponding certification). It permits only one sieving level. The mesh is easily interchanged and self-tensing.

It is available in 405 and 600 mm diameters.



### **Examples of application:**

Enamels, paints, ceramic powder, sand, seeds, chemicals, agriculture products, food, etc.





### **ADVANTAGES:**

- Mobile sieve, to work in the field.
- Highly versatile. Easy to operate.
- Interchangeable and self-tensing mesh.
- Adjustable height.
- Low maintenance.

### Tension (V) Intensity (A) Power (Kw) Speed (rpm) \*ØxH mm 50 Hz 60 Hz 50 Hz 60 Hz 50 Hz 60 Hz 50 Hz 60 Hz FTIS-0450 Circular 450 400 460 0,155 0,175 1500 1800 460 x 880 FTIS-0600 Circular 400 460 0,60 0,60 0,24 0,27 1500 1800 645 x 1000

### FTIC-0650 / FTIC-0800

The MINERVA centrifuge sieve shakers are used when high productivity is required, with two types of sieving being possible:

- Separation sieving of different granulometries in loose products with poor flow characteristics which are difficult to sieve due to their high content of humidity, grease or static.
- Safety sieving of the finished product in order to eliminate possible impurities.

Used in industrial sectors such as food, pharmaceuticals, chemicals etc., they are suitable for a wide range of products such as flours, sugars, ceramics, varnishes, pigments, plastics, etc. They can separate granulometries from 75 microns to 5 mm.

The product to be sieved is conveyed to the interior of the cylinder via the feed screw and is distributed over the whole sieve surface via centrifugal paddles. The velocity of the paddles creates a centrifugal force which pushes the product through the mesh. The largest particles or those with impurities are thrown out. The velocity and the paddles can be adjusted to the type of product.

They can be made to order in AISI 304 or AISI 316 stainless steel. If the centrifugal sieve shaker is going to be used in areas classified as being potentially explosive environments then it is manufactured accordingly and provided with the ATEX certificate.



### **Examples of application:**

Flour, sugar, chocolates, ceramics, paints, pigments, plastics, cement, plaster, etc.





### **ADVANTAGES:**

- Sieve shaker for large production volume. High performance.
- Specially designed for humid, greasy or electrostatically charged products.
- Highly robust. It can work for 24
  hours
- Easy interchange of cylinder sieves.
- Fixed or portable.
- Silent.
- Easy to clean. Low maintenance.

MODEL	Sieve Shaker type	Sieve Diameter (mm)	Tension	Intensity (A)	Power (Kw)	Speed (rpm)	Weight (Kg)	Dimensions LxWxH mm
FTIC-0650	Centrifugal	2 x Ø 220 x 294	380/400 V 50 Hz	3,21	1,4	750	205	1875 x 500 x 840
FTIC-0800	Centrifugal	2 x Ø 310 x 370	380/400 V 50 Hz	5,83	2,2	900	245	2200 x 600 x 945

Filtra CE

# **ATENEA**

# Accessories and spare parts

### FTMV-16 / 25 / 75 / 250 / 500

The ATENEA industrial V-shaped mixer is used in order to achieve perfect homogenization of solid with solid (in any proportion) and solid with liquid (10% maximum) products in powder or granulated form.

It has a V-shaped body with two loading inlets and one outlet, with a butterfly valve to facilitate the unloading of the mixture. Due to the V-shape axial currents are produced when the body is turned through 360°, and these currents separate and unite the product to be mixed. This, along with the radial mixing action, produces an extremely homogeneous result (in the proportion of 1 part in 10.000).

It is possible to add an intensifying bar within the interior which breaks down lumps and reduces the mixing time. Using this intensifying bar it is possible to add liquids in spray form.

The surface finish is mirror-polished with various grades of polish as required. The mixing time varies between 3 and 15 minutes, depending on the complexity of the mixture. It is available in various models for capacities of 16, 25, 75, 250, 300 or 500 litres of product.

### **OPTIONAL:**


Filtra

- Intensifier (shaker) / humidifier bar set • External electrical panel.
- Safety protective barriers.
- Finishes and qualities.





- Precision mixing. Maximum product homogenization.
- High quality materials and
- Free of interior angles or corners; there is no build up of
- · Easy to maintain.

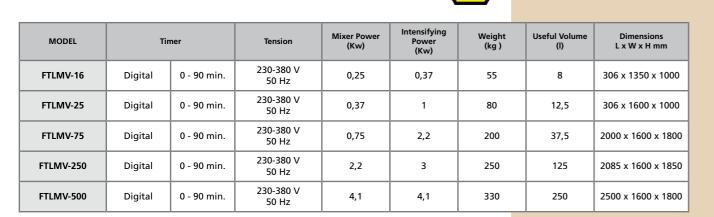

### **Examples of application:**

Flours and other food products in powder or grain, pigments, colourants, glass fibre, plastic in powder or granules, pesticides, herbicides, fertilizers, polymerized products, pharmaceutical products, chemical products, etc.





### **ADVANTAGES:**

- finishes.
- product residue.
- Silent.





### For Centrifugal Sieve Shaker:

- Mesh cleaning rings (4) - Quick clamp system (1)

- Sleeves: nylon, silicone (8)

- Silicone gaskets

- Remaking sieves (5)

- Support bench

- Silicone Lids (3)

- Springs (2)

- Wheels (7)

- Control (10)

- Cylindrical Sieves: St. steel or nylon (6)

- Sieve: 550, 800, 1200, 1500, 2000 mm. (5)

- Support bench
- Control (10)

### For "V" Mixer:


- Safety protective barriers: cage, cells, etc. (9)
- Support bench
- Control (10)
- Inox Lids (12)
- Intensifier (shaker) / humidifier bar set (11)









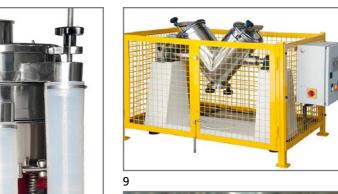

















FILTRA VIBRACIÓN, S.L.
Pol. Ind. Les Guixeres · C. Bronce, 1-3 · 08915 Badalona (Barcelona) · España / Spain
Tel. / Ph.: (+34) 933 956 739 · Fax: (+34) 902 190 070

filtra@filtra.com www.filtra.com